


Prehension trial progress and synchronization information

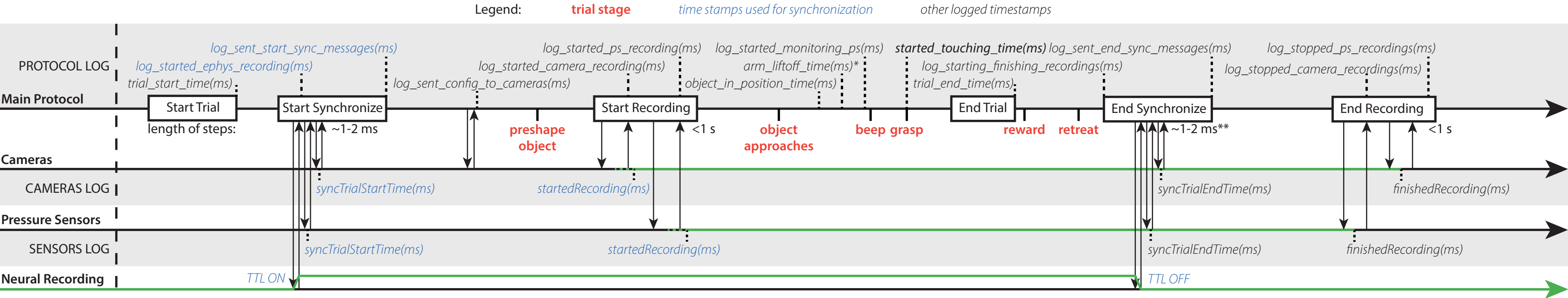

To align cameras or sensors to neural recording TTL ON pulse, add $startedRecording(ms) - syncTrialStartTime(ms)$ to their internal time of the trial start.

* Signal to start trial as shown happens after arm liftoff timestamp, however the arm liftoff is not actually tracked at this time. Rather, this timestamp marks end of a process watching early grasp and initialization of the process that watches success.

** Actual end synchronization time can be much higher right now, closer to 200 ms because the EPHYS is controlled through a USB port on one of the computers that records from cameras.